

Terahertz generator and detector






- **Efficient THz generation using optical rectification of femtosecond pump pulses**
- **Efficient THz generation using nonlinear optical difference frequency generation**
- **Optimized for pump wavelengths of 1.2–1.6 μm ; 0.7–0.8 μm available upon request**
- **Efficient electro-optic THz detector**

Specifications

Aperture	2, 3, 4 or 5 mm; others upon request
Damage Threshold	150 GW/cm ² @150 fs pulse length
Photon conversion efficiency	2·10 ⁻⁴ /MW-peak-power



THz Frequency Ranges for Generator Materials

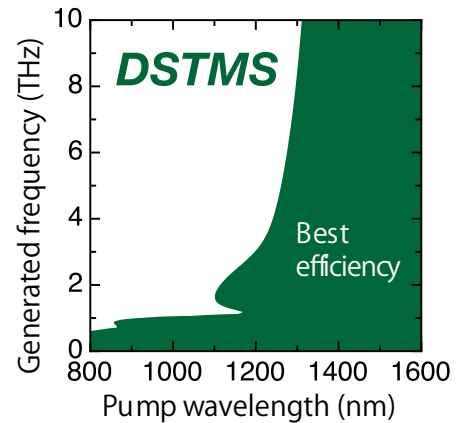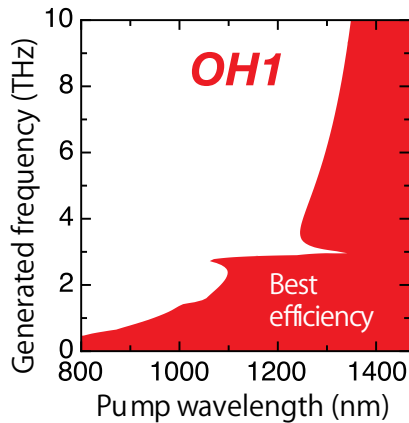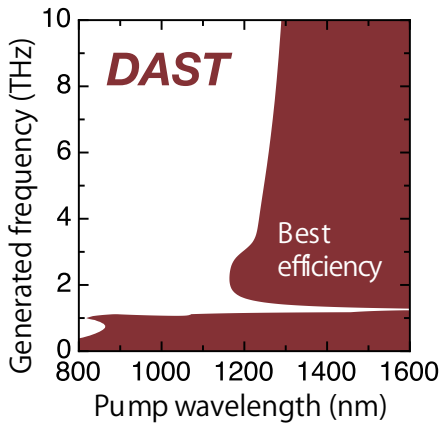


Spectral Bandwidth (measured with Rainbow Photonics instruments)

Source/Detector: 0.45 mm DSTMS
 $\lambda = 1560$ nm

Pump Pulse length: 65 fs
 Energy/Pulse: 1.8 nJ, Average Power: 180 mW

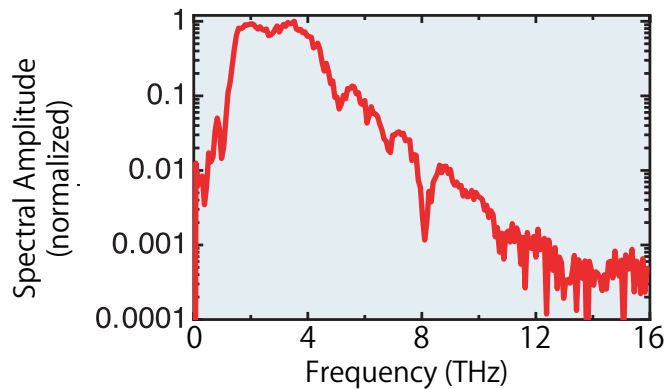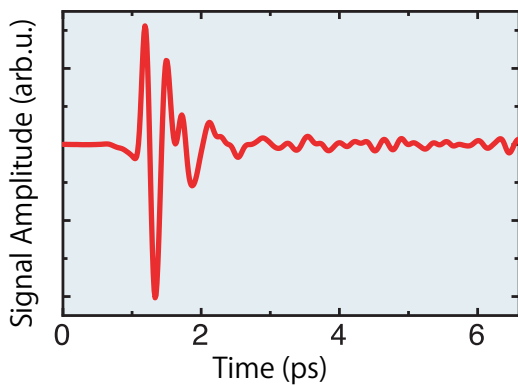

THz Frequency Range for Different Pump-Pulse Lengths

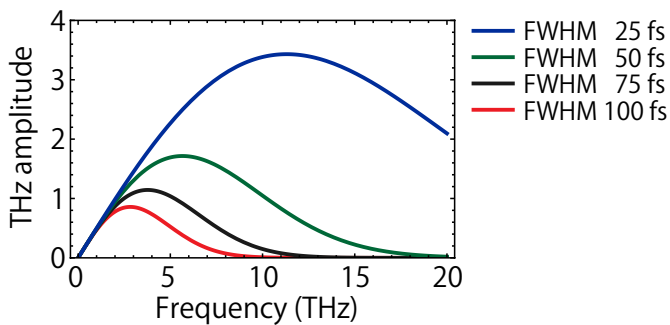

References

- A. Schneider et al, Appl. Phys. Lett. 84, 2229 (2004).
- A. Schneider et al, J. Opt. Soc. Am B 23, 1822 (2006).
- F. Brunner et al, Opt. Express 16, 16496 (2008).
- M. Stillhart et al, J. Opt. Soc. Am B 25, 1914 (2008).

More information available upon request.

Rainbow Photonics AG

Farbhofstrasse 21
 CH-8048 Zürich

Phone: ++41 44 419 05 05
 Fax: ++41 44 419 05 06
 E-mail: info@rainbowphotonics.com
 Web: www.rainbowphotonics.com



PNEUM Co., Ltd.

5-15-3 Minamikoshigaya, Koshigaya-shi,
 Saitama-ken, 343-0845, Japan

TEL: 81-48-985-2720
 FAX: 81-48-985-2721
 info@pneum.co.jp 1502